

Danis Procedure Pack

SX-ELLA Stent Danis



Danis



Removable stent system for treatment of bleeding esophageal varices

Easy implantation without the need for endoscopy / X-Ray control

Secure hemostasis through direct compression of varices

Oral intake maintained from time of implantation



3D animation



Basic description

Procedure Pack with SX-ELLA Stent Danis

Self-expandable nitinol stent preloaded in a ready-to-use delivery system. Complete procedure pack including all items required for implantation.

Indications

SX-ELLA Stent Danis is intended for stopping the acute bleeding from esophageal varices. The stent implantation may be an option in refractory esophageal variceal bleeding as an alternative method to the early TIPS or the balloon tamponade.

Contraindications

- ☑ Esophageal strictures
- ☑ Recent radiotherapy of the esophagus
- ☑ Tumors of the upper respiratory tract, esophagus or stomach
- ☑ Upper GIT bleeding from other causes (e.g. perforation)

Features / Benefits

- ☑ Variable pitches in stent weaving conform to esophageal peristalsis - reduced migration rates.
- ☑ Standardized compression of the varices results in effective hemostasis.
- ☑ Atraumatic stent ends.
- ☑ Radiopaque markers at both stent ends and at the midpoint - accurate stent positioning.
- ☑ Retrieval loops made of medical grade alloy at both stent ends - high mechanical durability and acid resistance.
- ☑ Easy endoscopic stent removal after 7 days.

Benefits of the method

- ☑ Readily implantable without endoscopy / X-Ray control, also in emergency situations.
- ☑ Patient can continue an oral diet after stent implantation.

Additional description


Danis Procedure Pack (019-08S-25-135) consists of a plastic carry case containing all items required for stent implantation:

- ☑ Delivery system preloaded with stent, guide wire, plastic syringes, mouthpiece, kidney dish, slobber cloth, waste bag with clip, sterile gloves.

Danis Procedure Pack Basic (019-08S-25-135-B) consists of:




- ☑ Delivery system preloaded with stent, guide wire, plastic syringe 50 ml.

The delivery system 9.4 mm (28 F) / 6.6 mm (20 F) is available in a standard length of 60 cm.

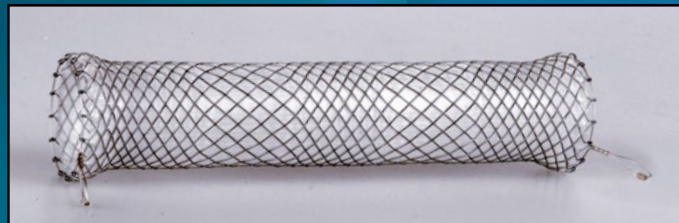

Available sizes

REF No.	SX-ELLA Stent Danis			Delivery system	
	Stent flares diameter [mm]	Stent body diameter [mm]	Nominal length [mm]	Active length [cm]	Outer diameter [F]
019-08S-25-135	30	25	135	60	28 / 20
019-08S-25-135-B	30	25	135	60	28 / 20

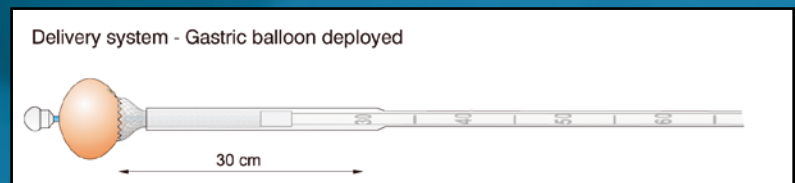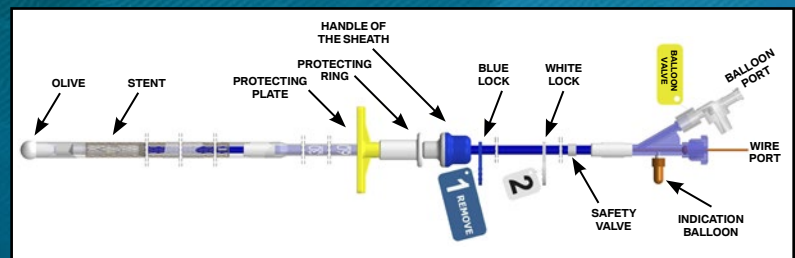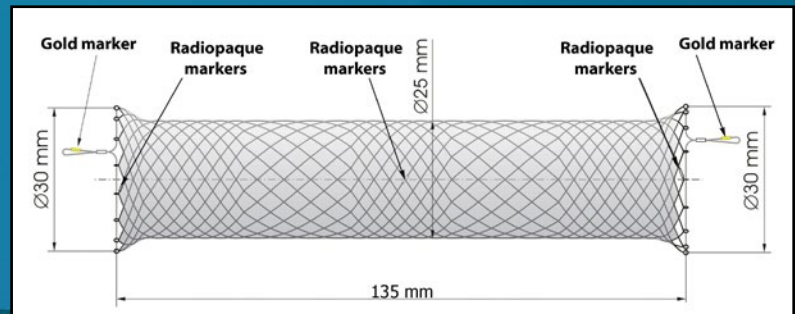

Distributor:

Manufacturer:

ELLA-CS, s.r.o.

Milady Horakove 504/45, Trebes
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655
E-mail: info@ellacs.eu URL: www.ellacs.eu